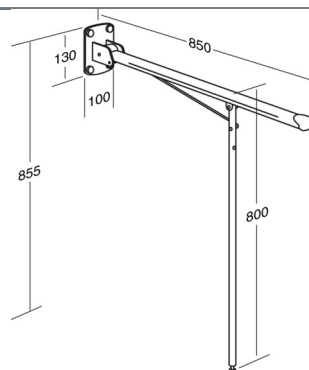




PRESSALIT®

| | |
|----------------------------|---|
| Product | VALUE support arm, fixed height |
| Item number | R1102000 |
| Colour | White |
| Product description | <p>VVS:615722380 NRF:6191677 RSK:7972411 EAN:5708590192345 HMI: 29695 NBS:N13:436 DB: 2109684 FDA:ILS/D135866</p> <p>VALUE support arm with leg, fixed, foldable, projection 850 mm, height 800 mm. The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p> <p>VALUE support arm with leg, fixed, foldable, projection 850 mm, height 800 mm. The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends. Materials: Rilsan-coated steel, pressure-cast aluminium, plastic. Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p> |
| Adjustability | Fixed mounting |
| Product capacity | <p>Maximum load: 135 kg. When the support arm is loaded with 135 kg, each wall fastener is loaded with 81 kg.</p> |
| Materials | The support arm is manufactured from steel with a Rilsan coating, which is comfortable to the touch, easy to clean, and protects the metal against corrosion. The end caps are manufactured from nylon (PA). |
| Ergonomy | <p>Provides firm support when you are sitting, standing, or transferring to and from a toilet.</p> <p>Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.</p> |





PRESSALIT®

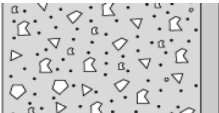
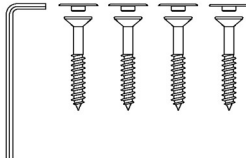

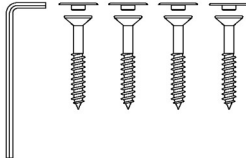
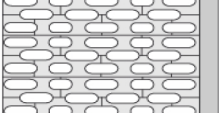
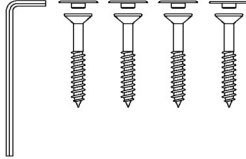
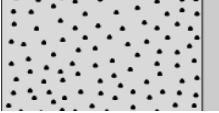
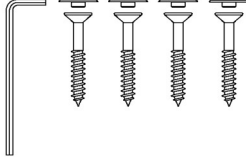
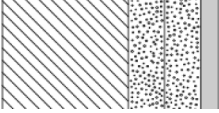
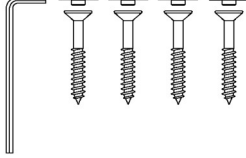
| | |
|----------------------------------|--|
| Product | VALUE support arm, fixed height |
| Item number | R1102000 |
| Features | <p>The support arm can fold up when not in use or to make more space.</p> <p>With leg for extra support.</p> |
| Warranty | Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date. |
| Cleaning | <p>VALUE support arms.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Klorilli</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> |
| Declaration of conformity | CE MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices. |



PRESSALIT®

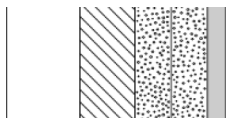
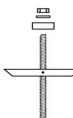

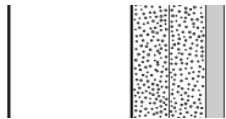
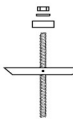

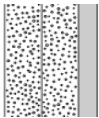

| | |
|--------------------|---------------------------------|
| Product | VALUE support arm, fixed height |
| Item number | R1102000 |

Installation in different wall types

| | | |
|---|--|---|
|  | <p>Concrete >B25</p> <p>V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> |  |
|  | <p>Solid Brick >= Mz12</p> <p>V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> |  |
|  | <p>Solid Sand Lime Brick >= Hlz12</p> <p>V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> |  |
|  | <p>Aerated Concrete P4</p> <p>V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> |  |
|  | <p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8593 - 4 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> |  |



PRESSALIT®

| Product Item number | VALUE support arm, fixed height R1102000 | | |
|---|--|---|---|
|  | Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar. |  |  |
|  | Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar. |  |  |
|  | Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls. |  | |

Bolts for wall fixing are included. It is the responsibility of the installer to ensure that suitable plugs are used. Please see mounting instruction.

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®